


# MALDI Imaging Mass Spectrometry for the study of cardiovascular pathology

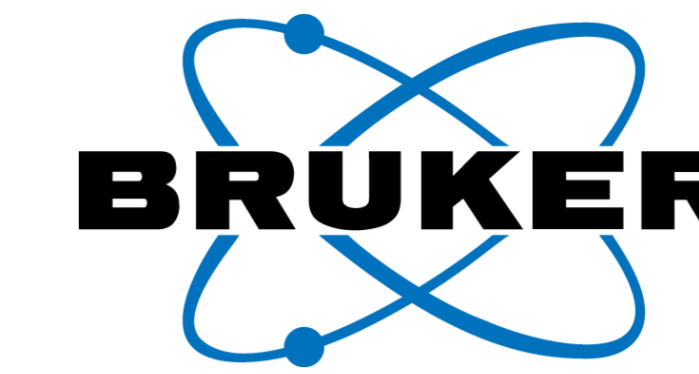

Takashi Nirasawa<sup>1</sup>, Megumi Terada<sup>2,3</sup>, Hiroko Namba<sup>2</sup>, Nobuto Kakuda<sup>2</sup>, Patrick Bruneval<sup>3</sup>, Hatsue Ishibashi-Ueda<sup>4</sup> and Masaya Ikegawa<sup>1</sup>

1. Bruker Japan, K.K., Yokohama, Japan 2. Doshisha University, Faculty of Life and Medical Systems, Kyoto, Japan 3. Georges-Pompidou European Hospital, Anatomy-Pathology, Paris, France 4. National Cerebral and Cardiovascular Research Center, Osaka, Japan

ASMS2018 TP276

## Overview

MALDI-IMS, Cardiovascular disease, Endomyocardial biopsy, FFPE

## Introduction

Endomyocardial biopsy (EMB) is performed on a regular basis following each transplant for the diagnosis of cardiac allograft rejection (CAR) based on the criteria established by International Society for Heart and Lung Transportation (ISHLT). Here we have applied a comprehensive Matrix-assisted laser desorption/ionization (MALDI) Mass spectrometry (IMS) approach to formalin-fixed paraffin-embedded (FFPE) tissue samples that have been acquired from the EMB conducted after the transplantations at Hôpital Européen Georges-Pompidou in Paris, France.

## Methods


**Tissue Preparation:** FFPE tissue samples were transferred from Georges - Pompidou European Hospital, (GPEH) Paris, France. Tissue sections were prepared from FFPE tissue samples according to modified procedures as follows: 1) Dewaxing 2) On-tissue digestion using trypsin as a protease 3) Deposition of MALDI matrix with Image Prep and HTX TM-sprayer with or without manual sprayer for better ionization.

**Immunohistochemistry:** Anti TTR antibody and anti AL antibody (Santa Cruze Biotechnology, Inc.) were used for IHC. Amyloid staining were performed with Sirius red.


**Mass Spectrometry:** Matrix coating: DHB in 50% MeOH, ImagePrep and or manually sprayed with airbrush; Standard trypsin deposition protocol; MALDI: Rapiflex; flexImaging 5.0 software. (Bruker Daltonics). Tissues were stained by H&E. PCA and pLSA analysis of mass spectra was performed with SCiLS Lab 2016B software. (Bruker Daltonics).

**Cardiac Amyloidosis** represents a growing number of diverse and devastating diseases which comprise cerebral forms such as Alzheimer's disease and cerebral amyloid angiopathy (Kakuda et al., 2017). Peripheral disorders including senile systemic, dialysis-related, primary (AL) and secondary (AA) amyloidosis. These diseases are characterized by the aggregation of normally soluble, natural proteins or their fragments into highly ordered amyloid fibrils that deposit in target organs leading to dysfunction.


## Results and Discussions

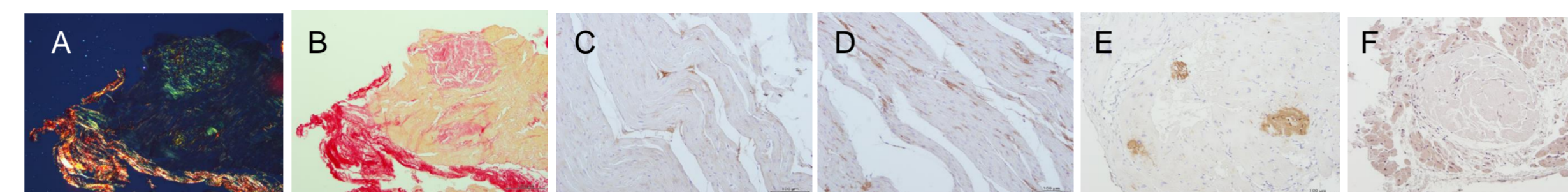


**Figure.1 Histochemistry of cardiac amyloidosis.** A,B: Sirius Red staining of cardiac amyloidosis. C,D: Cardiac fibrosis, E: TTR amyloidosis, F: AL amyloidosis. C,E: anti TTR antibody, D,F: anti AL antibody.

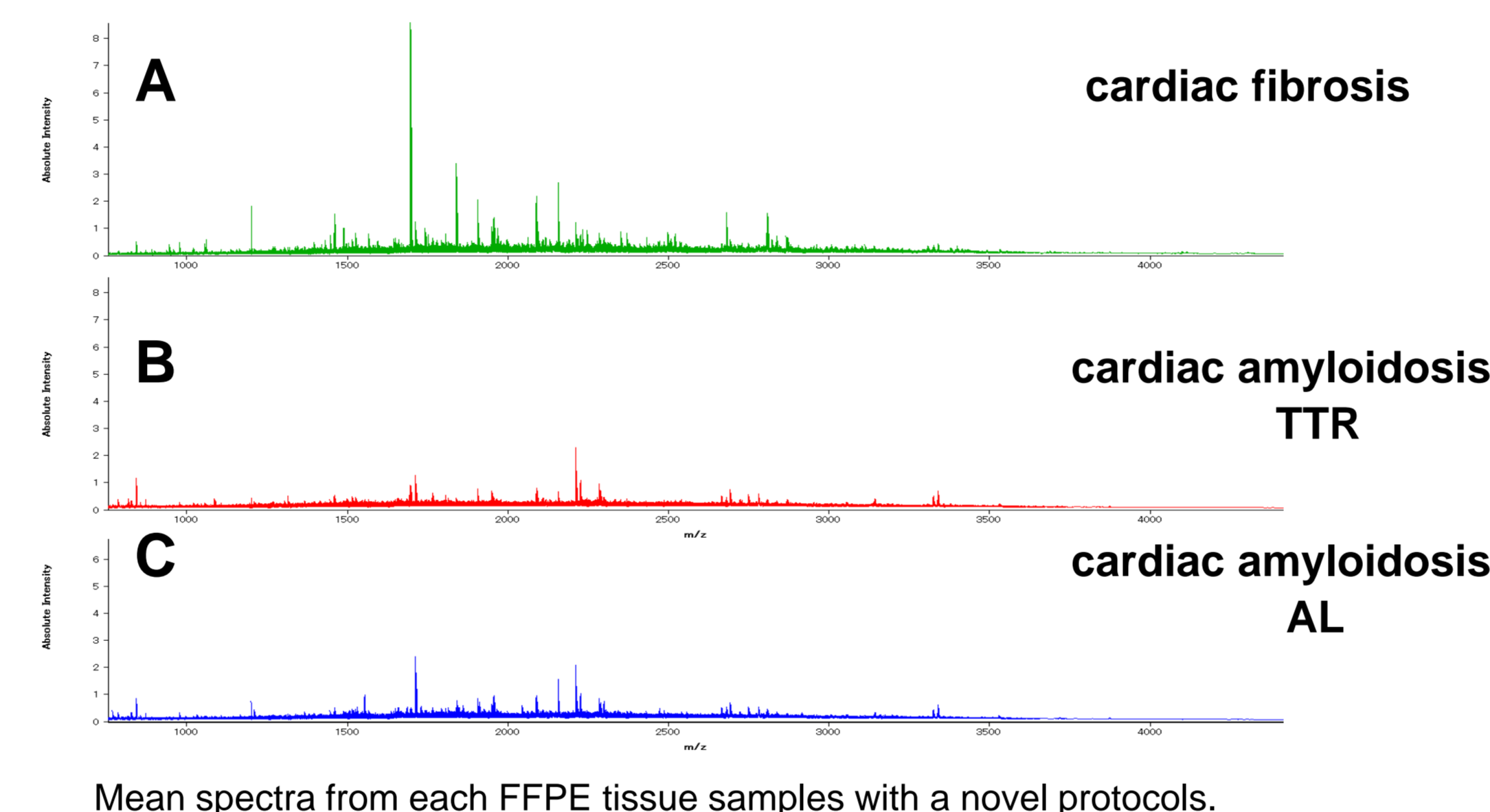

**Figure.2 Mass spectrometry obtained from FFPE samples.** Mean spectra from EMB samples from patients with cardiac amyloidosis and fibrosis was shown. A: Cardiac fibrosis, B: TTR amyloidosis, C: AL amyloidosis.

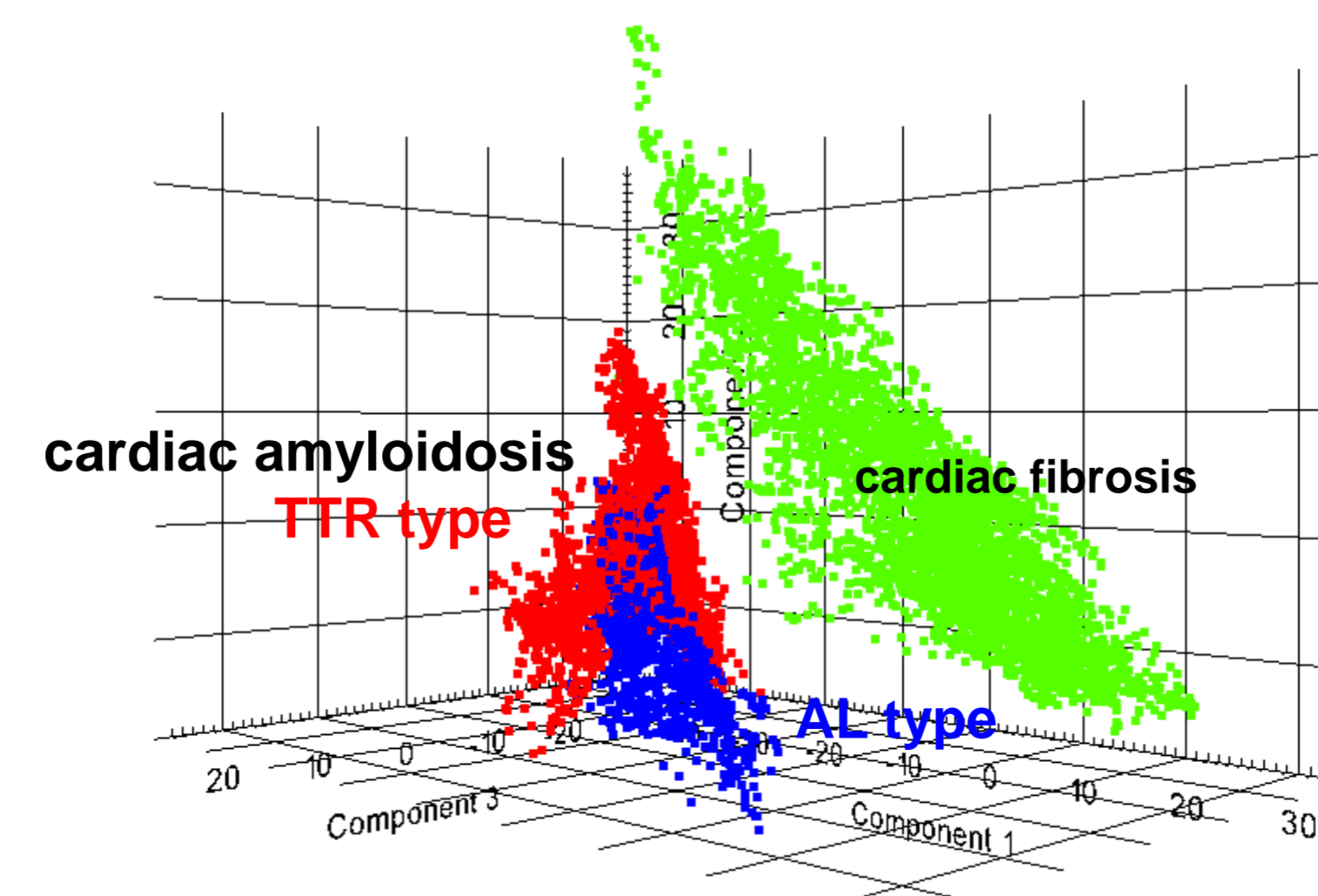

**Figure.3 PCA clearly differentiates cardiac fibrosis as well as subtype of cardiac amyloidosis.**

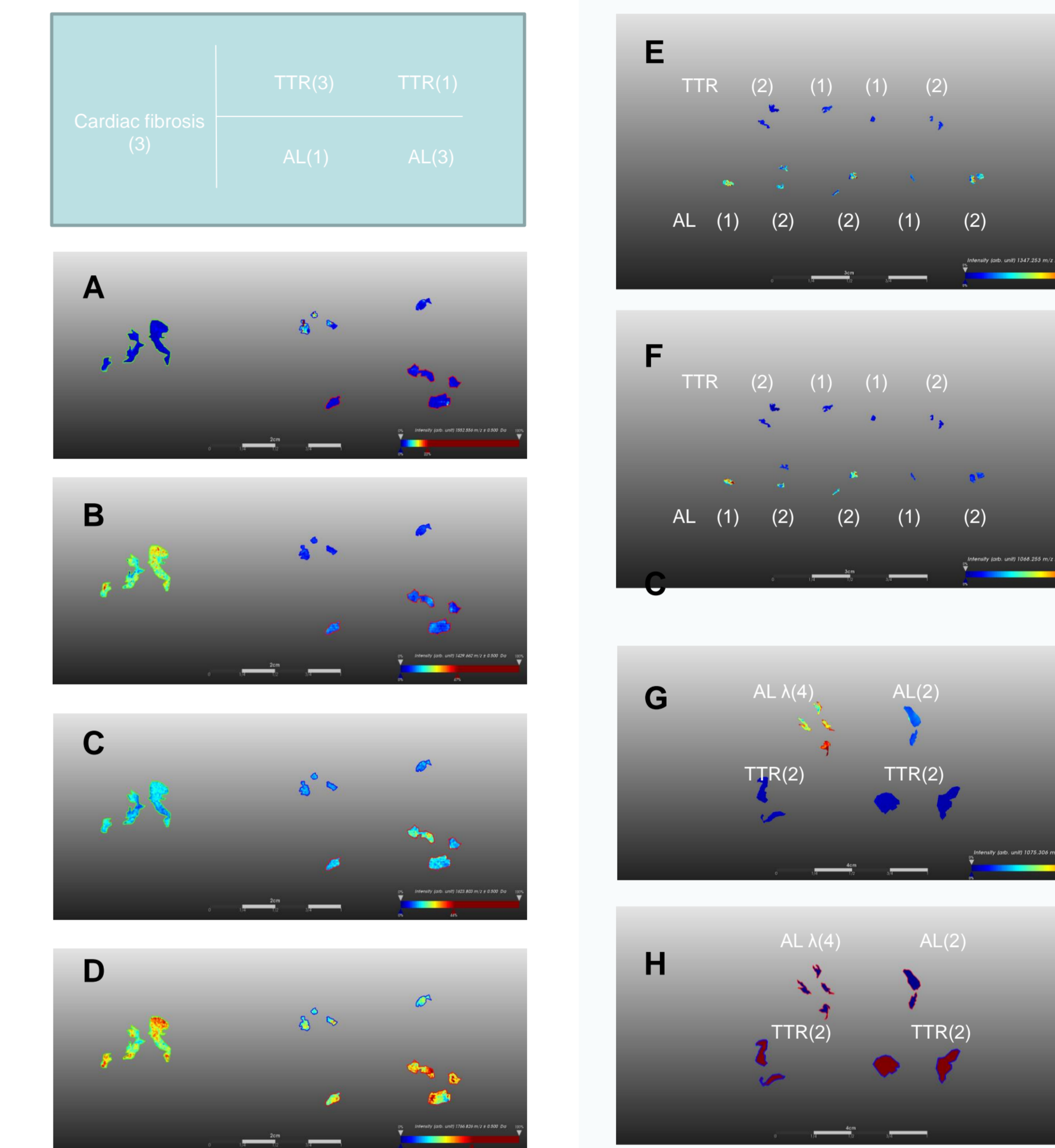

**Figure.4 MS ion image of each peptide in biopsied samples.** From A to D, A. TTR type. B. Fibrosis. C. AL type D. AL + Fibrosis dominant single peak distribution was shown. From E to H: matrix spraying with HTX TM-sprayer. G. m/z 1075.3 derived from YASWYQQK, from 48th to 56th positions of a tryptic digest of human immunoglobulin  $\lambda$ . Ref.3

## Summary

- MALDI-IMS protocol for FFPE samples of human cardiac biopsy was established
- Typing cardiac amyloidosis from FFPE samples with MALDI-IMS
- Future studies using cardiac biopsy samples from cardiac allograft, atherosclerosis and cardiomyopathy is ongoing

## Reference

- 1) Circulation. Lupy A. Et al., 2017
- 2) Acta Neuropatho Comm. Kakuda et al., 2017
- 3) Clinical Biochemistry. Nakanishi et al., 2013
- 4) Circulation Reserach Nour-Eldine W., 2018

## Acknowledgements

We thank Miss Yuki Kuzuhara and Yume Mukasa for experimental assist and scientific contributions for the establishment of FFPE protocol.