

## Application Note SC-XRD 505

# Single Crystal Diffraction

### Introduction

Single-crystal X-ray diffraction, commonly referred to as X-ray crystallography, is an analytical technique in which X-ray methods are employed to determine with certainty the actual arrangement of atoms within a crystalline specimen. The science of X-ray crystallography originated in 1912 with the discovery by Laue that crystals diffract X-rays. Since that time, single-crystal X-ray diffraction has developed into the most powerful method known for obtaining the atomic arrangement in the solid state.

X-ray crystallographic structure determination can be applied to a wide range of sizes of structures, from very small molecules and simple salts, to complex minerals, synthetically prepared inorganic and organometallic complexes, natural products and biological macromolecules, such as proteins and even viruses.

The precise knowledge of the molecular geometry is important in nearly all fields of chemical and biological research. The three-dimensional atomic coordinates obtained from crystallographic studies, available in comprehensive crystallographic databases such as the Protein Data Bank and the Cambridge Structural Database, are often the starting point for molecular modeling, drug design and molecular orbital calculations. Indeed, many of the most significant advances in structural chemistry and structural biology are based upon results obtained from X-ray crystallographic analyses.

The results of X-ray crystallographic analyses are complementary to other commonly used solid-state techniques, such as X-ray powder diffraction, solid-state NMR, EPR, FT-IR, Raman spectroscopy, and neutron diffraction. Chemists also routinely use such techniques as nuclear magnetic resonance, infrared and ultraviolet spectroscopy, mass spectrometry, X-ray fluorescence, and elemental analysis for the identification and characterization of compounds prepared and isolated in their laboratories. The experimental data obtained from these techniques may, after suitable analysis and interpretation, yield important information concerning the composition and structure of the compound. However, such information is often incomplete, fragmentary and ambiguous. There are many classes of chemical compounds such as natural products, organometallic complexes, inorganic salts, metal cluster systems, and organic and inorganic reaction products for which the structure cannot be deduced even with all of the other methods combined. X-ray crystallography is uniquely capable of unambiguously determining the complete three-dimensional molecular structures (including the absolute stereochemistry) of chemical substances. Modern X-ray crystallographic data often permits routine location and refinement of solvent molecules and hydrogen atoms.

A significant reason for the increased use of X-ray crystallographic methods is the remarkable improvement in instrument performance achieved through the introduction of two-dimensional detectors, most recently CMOS systems. The improvements in X-ray crystallographic instrumentation, along with associated advances in computer hardware and software, have resulted in enormous progress in efficiency and productivity. Crystal screening and identification may now be performed in a matter of minutes, while analytical crystal structures may be completed in far less than one hour and publication-quality structures are typically done in less than 3 hours, depending upon the diffraction properties. Crystallographic calculations that once required overnight execution on a mainframe computer are now easily carried out on a notebook computer in minutes. A single instrument may now be used to screen and identify many thousands of samples and to determine hundreds of structures per year, dramatically reducing the cost per structure.

X-ray crystallography, once regarded as an expensive and time-consuming technique used only by specialists, continues to gain new users in all branches of chemistry and biochemistry. Modern commercial instruments feature easy-to-use graphical user interfaces and automated crystallographic routines that allow routine crystal structure analyses to be carried out quickly and easily by users with only minimal training in X-ray crystallography. The high throughput and low operating costs of modern single-crystal instruments and their ease-of-use now make them suitable for use in synthetic laboratories as routine analytical tools.

The improvements in speed of data collection and data quality offered by new instrument technology have also extended the use of X-ray crystallography to a further level. Charge-density studies, which were previously limited to very small molecules due to long data collection times

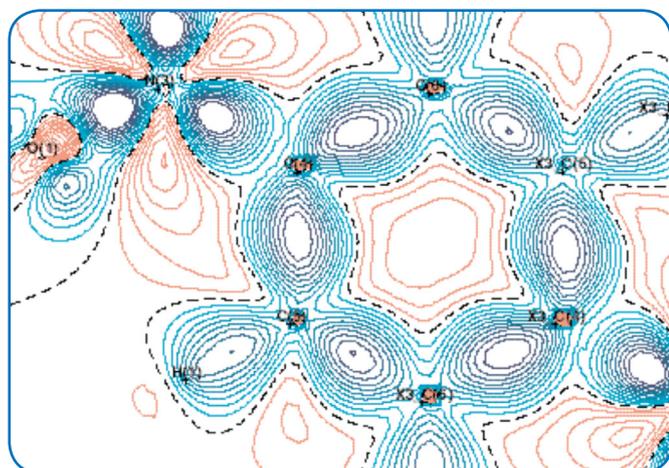


Figure 1: Electron density difference map for *p*-dinitrobenzene.

(often 3 to 6 months on conventional instruments), may now be carried out in a couple of days on much larger molecules. At this level, it becomes possible not only to determine with very high positional accuracy the three-dimensional positions of atoms in molecules and hence their bond lengths and angles, but also a measure of bond strength, atomic charges, dipole moments and the electrostatic potential. All of these properties are extremely important in predicting the chemical behavior of molecules (Figure 1).

The range of materials that may be studied by X-ray crystallography has recently been extended by the availability of higher-flux X-ray sources (e.g., microfocus sources, rotating anode generators, X-ray optics, liquid-metal sources, synchrotron beam lines) and more sensitive X-ray detectors (CMOS detectors), by the use of high-speed computers with large amounts of mass storage and the development of new algorithms for solving and refining large and problematic structures (e.g. twinned specimens, incommensurate structures). These advances, coupled with progress in the related fields of crystal growth, cryo-crystallography, and crystal-mounting techniques, now permit X-ray structure determinations to be carried out on very small specimens (minimum dimensions of a few microns), on materials with very large unit cells (e.g. axis length 400 Å) and on materials that are liquids at room temperature or that undergo solid-state phase changes.

## Fundamentals

A perfect crystalline solid (single crystal) is made up of a large number of identical molecules which are arranged in a precisely regular way repeated in all directions to give a highly ordered structure. The basic building block in a crystal is the unit cell. A crystal is made up of millions of identical unit cells arranged in a three-dimensional crystal lattice. Each crystalline substance has a unique set of **lattice constants** ( $a$ ,  $b$ ,  $c$ ,  $\alpha$ ,  $\beta$ ,  $\gamma$ ) which define the size and shape of the unit cell (Figure 2).

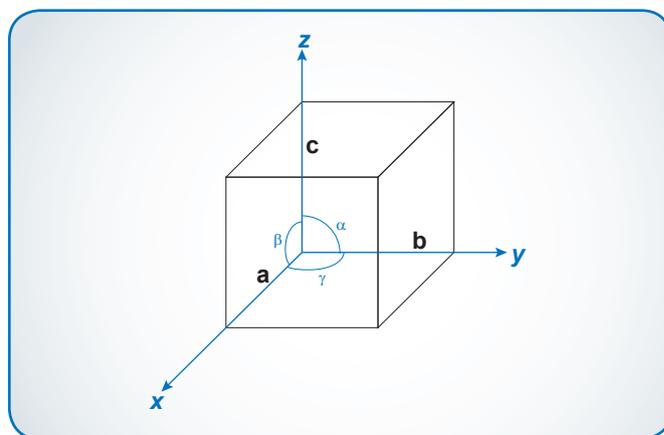


Figure 2: Unit-cell nomenclature: the reference axes  $x$ ,  $y$ ,  $z$  are right-handed. The length of the unit cell edge parallel to each reference axis is  $a$ ,  $b$ ,  $c$ , respectively. The inter-axial angles are  $\alpha$ ,  $\beta$ ,  $\gamma$ , respectively.

Based on the lattice constants and the symmetry that the unit cells possess, the substance may be classified into one of **seven crystal systems** (i.e., triclinic, monoclinic, orthorhombic, tetragonal, cubic, trigonal or hexagonal). Finally, each substance may be further classified as belonging to one of **230 three-dimensional space groups**. The lattice constants, crystal system, and space group are important physical constants for crystalline substances and may be used in conjunction with other physical measurements (e.g., density, conductivity, hardness) to explain the properties of solid-state materials. When a beam of parallel monochromatic X-rays of approximately one Angstrom wavelength strikes a single crystal, the crystal acts as a three-dimensional diffraction grating and produces an X-ray diffraction pattern (Figure 3).

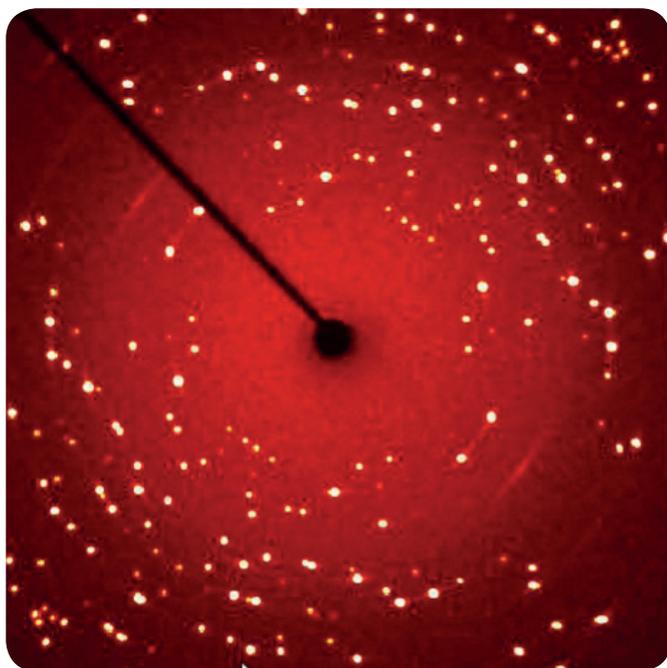


Figure 3: An example of an X-ray diffraction pattern (one-minute rotation image) produced by a randomly oriented single-crystal specimen on a 2D detector.

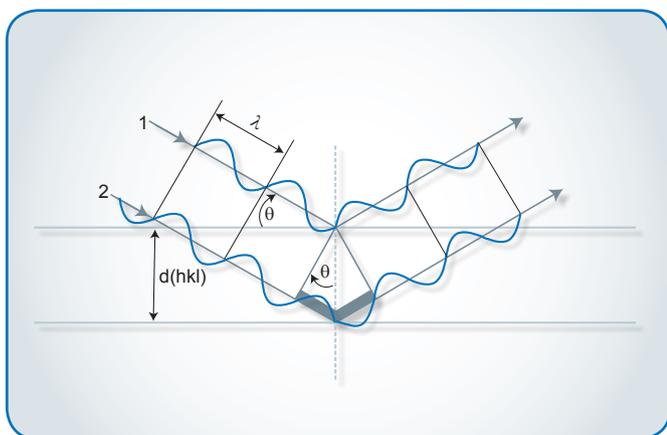


Figure 4: Diagrammatic representation of Bragg's Law showing the diffraction angle  $\theta$  and the interplanar spacing  $d(hkl)$ .

This diffraction consists of a three-dimensional array of reflections that satisfy the conditions of Bragg's law:

$$n\lambda = 2d \sin\theta$$

where  $n$  is a small integer giving the order of diffraction,  $\lambda$  is the wavelength of the incident X-rays,  $d$  is the distance between a set of parallel lattice planes, and  $\theta$  is the angle between the incident X-ray beam and the atomic lattice plane in the crystal dimensions of the unit cell in the crystal (Figure 4).

Each reflection may be assigned a set of indices (hkl) which indicate its location in the diffraction pattern or reciprocal space. The diffraction pattern in reciprocal space has a Fourier transform relationship to the electron density in the unit cell in real space.

Experimentally, the unit-cell parameters for a crystalline specimen may be determined from an analysis of the spatial arrangement of the reflections in its X-ray diffraction pattern. On modern instruments, the measurement of reflection positions and the subsequent calculation of unit cell parameters is an automatic process which is carried out in a few minutes as part of the specimen screening process. There are now several commercial databases (e.g., CSD) that contain unit cell, space group, chemical composition and bibliographic data on hundreds of thousands of compounds for which X-ray structure determinations have been carried out. Experimental unit cell data obtained in the crystal screening process and chemical composition information for the sample being studied may be used to search the databases. If an exact match is found, the compound has been identified and further single-crystal analysis may be unnecessary.

The analysis of the positional information for reflections, often called geometric data collection, may be used to provide direct information concerning the size and shape of the cell. In combination with density and elemental analysis measurements, it may also yield information concerning the cell contents; however, it does not reveal the actual locations of the atoms within the unit cell.

The determination of the arrangement of atoms in the unit cell requires a very detailed analysis of the relative intensities of all unique reflections in the diffraction pattern. The precise measurement of the relative intensities of the reflection is termed intensity data collection. The collection of accurate intensity data requires a highly stable X-ray source, a precise mechanical goniometer for sample positioning and a very efficient counting system. An intensity data set consists of several thousand reflections indexed by **h**, **k** and **l** for which an integrated intensity **I(hkl)** has been accurately measured.

The X-ray diffraction pattern consists of the superposition of scattered waves of varying amplitude and phase. Each diffraction maximum or reflection has associated with it a structure factor  $\mathbf{F}(\mathbf{hkl})$  measured relative to the scattering by a single electron. The structure factor may be represented as a complex vector:

$$\mathbf{F}(\mathbf{hkl}) = \mathbf{A}(\mathbf{hkl}) + i\mathbf{B}(\mathbf{hkl})$$

where  $\mathbf{A}(\mathbf{hkl})$  and  $\mathbf{B}(\mathbf{hkl})$  are the real and imaginary components of  $\mathbf{F}(\mathbf{hkl})$  (Figure 5).

The magnitude or length of the vector  $|\mathbf{F}(\mathbf{hkl})|$  may then be represented as

$$\begin{aligned} |\mathbf{F}(\mathbf{hkl})| &= \{[\mathbf{A}(\mathbf{hkl}) + i\mathbf{B}(\mathbf{hkl})] \times [\mathbf{A}(\mathbf{hkl}) - i\mathbf{B}(\mathbf{hkl})]\}^{1/2} \\ &= [\mathbf{A}(\mathbf{hkl})^2 + \mathbf{B}(\mathbf{hkl})^2]^{1/2} \end{aligned}$$

Alternatively,  $\mathbf{F}(\mathbf{hkl})$  may be expressed as an exponential quantity:

$$\mathbf{F}(\mathbf{hkl}) = |\mathbf{F}(\mathbf{hkl})| \exp [i\alpha(\mathbf{hkl})]$$

where  $|\mathbf{F}(\mathbf{hkl})|$  is the amplitude of the scattered wave and  $\alpha(\mathbf{hkl})$  is its phase angle.

From Figure 5 it may be seen that

$$\mathbf{A}(\mathbf{hkl}) = |\mathbf{F}(\mathbf{hkl})| \cos \alpha(\mathbf{hkl})$$

$$\mathbf{B}(\mathbf{hkl}) = |\mathbf{F}(\mathbf{hkl})| \sin \alpha(\mathbf{hkl})$$

And, therefore

$$\tan \alpha(\mathbf{hkl}) = |\mathbf{B}(\mathbf{hkl})| / |\mathbf{A}(\mathbf{hkl})|$$

$|\mathbf{F}(\mathbf{hkl})|$  may be calculated directly from the measured intensity  $I(\mathbf{hkl})$  for a reflection, since

$$I(\mathbf{hkl}) = K|\mathbf{F}(\mathbf{hkl})|^2$$

where  $K$  is a constant. However, the phase angle  $\alpha(\mathbf{hkl})$  cannot be measured experimentally and must therefore be obtained indirectly through a variety of numerical techniques.

The central problem in the solution of a crystal structure is the assignment of phase angles to each reflection in the data set. The solution of the phase problem is considerably simplified for crystals that possess crystallographic centers of symmetry, since, to a first approximation, the imaginary components  $\mathbf{B}(\mathbf{hkl})$  are zero for centrosymmetric space groups and the phase angles are therefore restricted to values of  $0^\circ$  or  $180^\circ$ . A structure is considered solved when a set of phase angles has been found that allows the atoms

to be located and the calculated diffraction pattern to be matched to the experimental diffraction pattern.

Since the electron density in a crystal varies continuously and periodically in three-dimensional space, the electron density  $\rho(\mathbf{xyz})$  at a point with fractional coordinates  $\mathbf{x}, \mathbf{y}, \mathbf{z}$  in a unit cell of volume  $V$  may be expressed as a three-dimensional Fourier series:

$$\rho(x, y, z) = V^{-1} \sum_n \sum_k \sum_l |F_{hkl}| \cos 2\pi(hk + ky + lz - \alpha'_{hkl})$$

If both the amplitude  $|\mathbf{F}(\mathbf{hkl})|$  and the phase  $\alpha(\mathbf{hkl})$  of each reflection are known, the electron density within the unit cell of the crystal can be calculated directly. On the other hand, if the positions of the atoms in the unit cell are known, both the structure factor and the phase for each reflection may be calculated from the structure factor equation:

$$F_{hkl} = \sum_j^N f_j e^{2\pi i(hx_j + ky_j + lz_j)}$$

where  $f_j$  is the atomic scattering factor for the atom  $j$  and  $x_j, y_j, z_j$  are its fractional coordinates. In an actual structure determination both forms of the Fourier transform equations are utilized to arrive at a model structure from which the observed diffraction pattern can be reproduced.

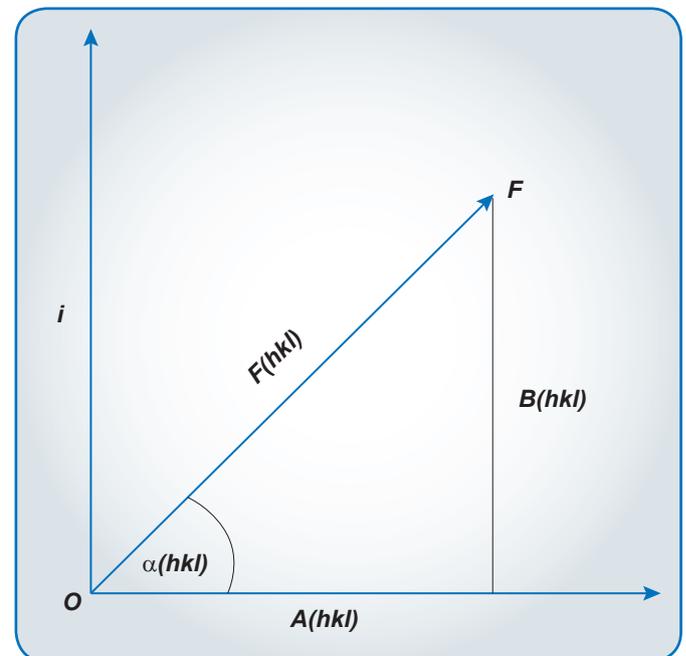


Figure 5: Structure factor  $\mathbf{F}(\mathbf{hkl})$  plotted on an Argand diagram.  $\alpha(\mathbf{hkl})$  is the phase angle and the amplitude is represented by length  $OF$ .

## Instrumentation

The basic hardware components of a typical automated single-crystal X-ray diffractometer system include:

- An X-ray source consisting of a high-stability X-ray generator, a copper or molybdenum or silver target X-ray tube, a tube shield with associated shutters, attenuators and safety interlocks, a monochromator or X-ray mirror system, and an incident-beam collimator.
- A three- or four-circle goniometer system that allows the specimen to be precisely oriented in any position while remaining in the X-ray beam.
- A video camera for aligning the specimen and indexing crystal faces.
- A CMOS-based two-dimensional X-ray detector system.
- A low-temperature attachment for cooling the specimen during data collection.
- A microprocessor-based interface module that receives commands from a host computer and carries out all real-time instrument control functions to drive goniometer motors, monitor the detector system, open and close the shutter and monitor collision sensors and safety interlocks.
- A host computer with a large hard disk mass storage device, a video monitor and keyboard, and diffractometer control programs to control the data collection strategy and to send commands to the microprocessor.

Structure determination calculations may be carried out on the computer used for data collection or they may be performed on a second computer linked to the diffractometer system. A large variety of hardware configurations are available, depending upon the requirements of the individual laboratory.

The most critical mechanical component in an X-ray diffractometer system is the goniometer assembly, which must be capable of keeping the specimen centered in the incident X-ray beam while at the same time changing its orientation in order to collect many thousands of frames of data in reciprocal space. The most commonly used type of goniometer is illustrated in Figure 6. The most important and most expensive component of any modern single-crystal diffractometer system is its detector system. Most new instruments purchased since 2011 use a detector system based upon 2D CMOS technology (Figure 7).

CMOS technology provides many advantageous features and is quickly developing. The PHOTON 100 CMOS detector (Figure 8) employs a large detection area of 100 x 100 mm, without need for a fiber-optics taper to increase the detection size. This improves the x-ray efficiency and increases the total number of X-ray reflections collected per image.

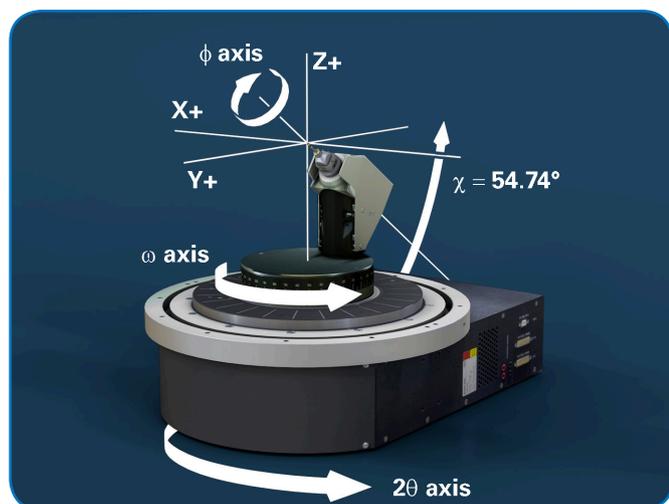


Figure 6: Diagram of axes for a three-circle goniometer ( $\chi$  is fixed at 54.74°).



Figure 7: Commercial CMOS-based single-crystal X-ray detector



Figure 8: A commercial CMOS-based single-crystal X-ray diffractometer system

## Experimental Procedure

The first step in a crystal structure analysis is concerned with the selection and mounting of a suitable specimen. Ideally, a crystal whose structure is to be determined must be a single crystal of 0.1 mm to 0.5 mm size, not cracked and not twinned. Micron-sized crystals are not uncommon with the higher-intensity microfocus sources. The techniques required to obtain such crystals vary considerably depending on the types of compounds to be analyzed. Stable crystals of typical organic, organometallic or coordination complexes can usually be grown by slow recrystallization from common solvents. Other types of compounds may require the use of sublimation, zone refinement, solvent-diffusion, low-temperature and/or inert-atmosphere techniques in order to isolate suitable specimens.

- Once a suitable specimen has been selected, it is glued or otherwise securely attached to a goniometer head (sample holder) in an arbitrary orientation.
- The goniometer head is then placed on the base of the goniometer assembly and the crystal is optically aligned in the center of the incident X-ray beam using a video camera or microscope. The orthogonal X-, Y-, and Z-translations on the goniometer head are adjusted until the specimen is centered for all crystal orientations.
- A preliminary rotational image is then collected with the 2D detector to screen the specimen for analysis and to select suitable parameter values for subsequent steps.
- In order to determine the unit cell, a preliminary set of frames is measured using an automatic routine. For example, three sets of frames are collected in different parts of reciprocal space.
- These frames are then processed to locate spots on individual frames and to then determine the centers of reflections.
- An auto-indexing routine selects the appropriate reduced primitive unit cell and calculates the corresponding orientation matrix and lattice constants.
- This preliminary unit cell is then refined using a non-linear least-squares algorithm and converted automatically to the appropriate crystal system and Bravais lattice. This new cell is refined by the non-linear least-squares algorithm to yield an accurate orientation matrix which may be used to index crystal faces and to carry out integration calculations after intensity data collection.

After the above geometric data collection steps have been completed and an accurate orientation matrix has been calculated, intensity data collection is carried out. Typically, a sphere or hemisphere of data is collected using a narrow-frame scan method in which several sets of frames (runs) are collected by scanning in  $0.1^\circ$  to  $0.3^\circ$  increments in the  $\omega$  and/or  $\phi$  angle, while keeping all other instrument angles constant. There are options to limit data to a unique set of reflections, thereby reducing data collection times for high-symmetry crystal systems. A complete data set may require anywhere from minutes to overnight depending on the size of the specimen and its diffracting power. Typical exposure times of 10 to 30 seconds per frame for a hemisphere of data require only a couple of hours of total data collection time, and may be as quick as 10 minutes for 1 second exposures in shutterless mode on the CMOS detectors.

When the complete set of frames has been collected for a given specimen, the entire data set must be processed to obtain accurate integrated intensities for individual reflections. This process includes corrections for instrumental factors, polarization effects, X-ray absorption and possibly crystal decomposition. The integration process reduces the raw frame data, which require from 500 to 2000 megabytes of disk storage, to a small set of individual integrated intensities for each reflection. The final unit-cell constants are calculated from the centroids of many thousands of reflections selected from the entire data set and typically have relative errors of less than 1/10,000.

Once the structure amplitudes are known, the phase problem must be solved to find a self-consistent set of phases that can be combined with the structure factor amplitudes to obtain the electron density and thereby determine the structure of the crystal. A number of crystallographic techniques exist for obtaining the phases of diffracted waves; the most widely utilized approaches to the solution of phase problem involve the use of either vector methods based on  $|F(hkl)|^2$  or direct or statistical methods. Typically, the solution to the structure yields only a partial or approximate model, which must be improved by successive applications of Fourier-transform methods before the complete structure has been determined.

After the entire molecular structure has been determined, the approximate positions of the atoms are refined by non-linear least-squares techniques to give the best fit between the calculated and observed intensity data for the specimen. The refinement process yields very accurate values for atomic positions from which bond lengths, bond angles and other structural parameters may be calculated. Finally, upon completion of the X-ray diffraction analysis, the structure of the molecule or crystal lattice may be displayed or plotted (Figures 9 and 10).

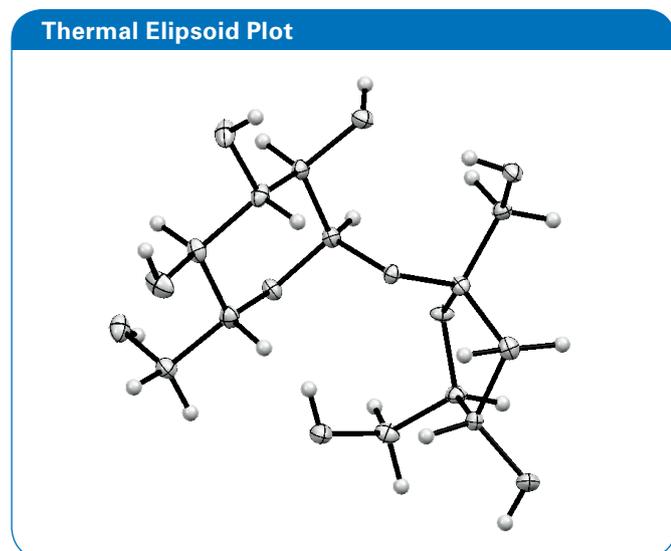


Figure 9: Thermal ellipsoid plot of the final structure of an organic compound ( $C_{12}H_{22}O_{11}$ ).

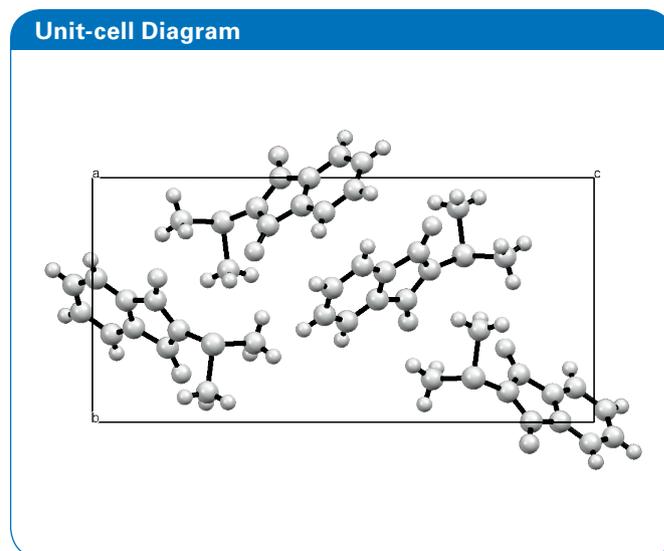


Figure 10: Unit-cell diagram showing the arrangement of molecules within the cell.

## Bibliography

**Crystal Structure Analysis.** Principles and Practice.  
(Second Edition)  
Clegg W, Blake A J, Cole J M, Evans J S O, Main P,  
Parsons S and D Watkin D J  
Oxford University Press, 2009

**Crystal Structure Determination.** (Second Edition)  
Massa W and Gould R O  
Springer, 2004

**Crystal Structure Analysis.** A Primer. (Third Edition)  
Glusker J P and Trueblood K N  
Oxford University Press, 2010

**Fundamentals of Crystallography.** (Third Edition)  
Giacovazzo C, Monaco H L, Viterbo, D, Scordari F, Gilli G,  
Zanotti G and Catti M  
Oxford University Press, 2011

**Crystal Structure Refinement**  
Müller, P., Herbst-Irmer, R., Spek, A.L., Schneider, T.R.  
and Sawaga, M.R.  
Oxford University Press, 2006

## Author

Charles Campana, Bruker AXS Inc.

---

 **Bruker AXS GmbH**

Karlsruhe · Germany  
Phone +49 721 50997-0  
Fax +49 721 50997-5654  
info.baxs@bruker.com

[www.bruker.com](http://www.bruker.com)