




S1 TITAN


- Integrated camera & small spot

• S1 TITAN integrated camera and small spot options


The S1 TITAN can be equipped with an integrated camera to provide visualization and accurate positioning of the measurement spot. The small spot option provides a small measurement area for the isolation of small features to be tested. Thanks to the S1 TITAN's SharpBeam™ optimized geometry, the precision and accuracy of the measurement are the same as for the normal spot- there is no need to extend the measurement time to achieve the desired precision. These options are available for both the S1 TITAN^{LE} and S1 TITAN^{SP}.

- Small spot isolates specific sampling area
- Camera ensures accurate measurement positioning
- Save up to 5 images per assay
- Provides record of measurement spot
- Images easily import into reports
- Camera resolution: 640 x 480

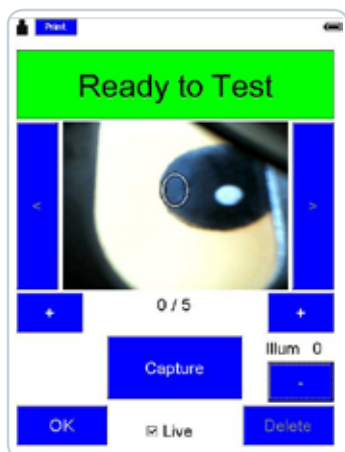



 Bruker-Elemental
 415 N. Quay
 Kennewick, WA 99336
 (This is editable, as is the logo)

ROHS_FP-35.pdz AssayTime: 7/9/2013 4:42:07 PM ElapsedTime: 21



3 of 5



1 of 5

Field Info

| | |
|--------|--------|
| Name | ID |
| Field1 | Field2 |

| Element Name | Min | PPM | Max | +/- [%] |
|--------------|---------|-------|-------------|---------|
| Cl | 0.000 | 32729 | 0.000 | 0.562 |
| Ca | 0.000 | 0 | 0.000 | 0.127 |
| Ti | 0.000 | 327 | 0.000 | 0.046 |
| Cr | 0.000 | 0 | 1000000.000 | 0.011 |
| Mn | 0.000 | 0 | 0.000 | 0.015 |
| Fe | 0.000 | 2 | 0.000 | 0.000 |
| Ni | 0.000 | 2 | 0.000 | 0.000 |
| Cu | 0.000 | 1 | 0.000 | 0.003 |
| Zn | 0.000 | 1 | 0.000 | 0.001 |
| As | 0.000 | 0 | 0.000 | 0.000 |
| Se | 0.000 | 0 | 0.000 | 0.000 |
| Br | 300.000 | 0 | 0.000 | 0.000 |
| Cd | 70.000 | 0 | 0.000 | 0.001 |
| Sn | 0.000 | 6 | 0.000 | 0.001 |
| Sb | 0.000 | 0 | 0.000 | 0.002 |
| Ba | 0.000 | 0 | 0.000 | 0.001 |
| Hg | 700.000 | 0 | 1300.000 | 0.01 |
| Pb | 700.000 | 0 | 1300.000 | 0.000 |

• Contact Us

www.bruker.com/s1titan

North America

Bruker
 Kennewick, WA · USA
 Tel. +1 (509) 783-9850
 hhsales@bruker-elemental.net

Europe / Middle East / Africa

Bruker
 Kalkar · Germany
 Tel. +49 2824 976 50-0
 hhsales@bruker-elemental.net

Asia / South America / Rest of World

Bruker
 Kennewick, WA · USA
 Tel. +1 (509) 783-9850
 hhsales@bruker-elemental.net